



# TriLux ICE

The Ultimate in Comfortable and Effective Hair Removal

At Aesthetics Consults, we are proud to present TriLux ICE, a revolutionary diode laser device for permanent hair reduction and skin rejuvenation. Designed with cutting-edge technology, TriLux ICE combines comfort, precision, and efficiency to deliver outstanding results for all skin tones and hair types



# Why Choose TriLux ICE?

TriLux ICE offers the perfect balance of advanced technology and user-friendly features, making it ideal for clinics aiming to provide superior hair removal treatments with maximum patient comfort.

## Key Features:

- **Three Wavelengths in One Handpiece:**

Targets different tissue depths and hair follicles for optimal results

- **Powerful Cooling System:**

Ensures pain-free treatments, even for sensitive skin

- **High-Speed Operation:**

Large spot size (12mm x 24mm) and up to 10Hz repetition rate for quick treatments

- **Safe for All Skin Tones:**

Customisable settings suitable for Fitzpatrick I–VI.

- **Durable and Reliable:**

Engineered for long-term use with minimal maintenance



# Comprehensive Treatment Capabilities

TriLux ICE is designed to address a variety of patient needs, offering more than just hair removal:

## PERMANENT HAIR REMOVAL

- Effective for fine, coarse, and light-coloured hair
- Treats all areas of the body, including face, underarms, back, chest, legs, and bikini line.

## SKIN REJUVENATION

- Improves overall skin texture and tone
- Boosts collagen production for smoother, firmer skin

## PIGMENTATION AND VASCULAR THERAPY

- Corrects pigmentation issues, including sunspots and melasma
- Reduces redness, rosacea, and small vascular lesions



# How TriLux ICE Works

TriLux ICE leverages three wavelengths and advanced diode laser technology for superior performance:

## Key Technologies:

### THREE-WAVELENGTH PRECISION:

- **755nm:** Ideal for lighter hair and superficial follicle targeting
- **808nm:** The gold-standard wavelength for most skin and hair types
- **1064nm:** Safe for darker skin tones with deeper penetration for coarse hair

### ICE COOLING TECHNOLOGY:

- Ensures constant skin cooling during treatments for maximum comfort.
- Prevents overheating and reduces the risk of skin irritation

### HIGH REPETITION RATE:

- Allows for faster treatment sessions without compromising results

TriLux ICE





# Benefits of Adding **TriLux ICE** to Your Clinic

## FOR YOUR CLINIC:

- **Expand Your Services:** Offer advanced hair removal and skin rejuvenation treatments.
- **Increase Revenue:** Treat more patients in less time with high-speed functionality
- **Reliable Performance:** Designed for long-term use with minimal downtime
- **Versatile Applications:** Safe and effective for all skin types, including sensitive and darker skin tones

## FOR YOUR PATIENTS:

- **Pain-Free Treatments:** ICE cooling technology ensures a comfortable experience
- **Quick Sessions:** Complete a full leg treatment in under 20 minutes
- **Visible Results:** Permanent hair reduction after 6–8 sessions
- **No Downtime:** Patients can return to daily activities immediately

## Hair Removal:

Noticeable reduction in hair growth after just a few sessions.




Before




After

## Skin Rejuvenation

Smoother, firmer skin with improved elasticity and tone

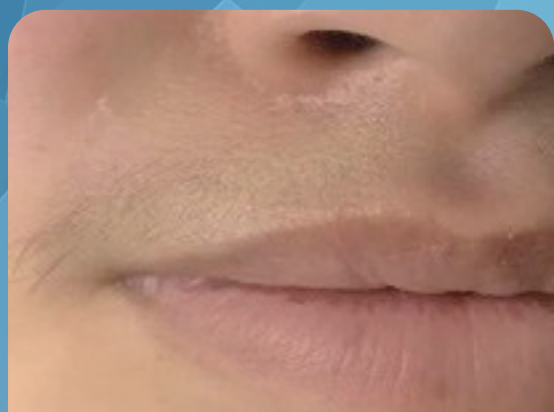

Before

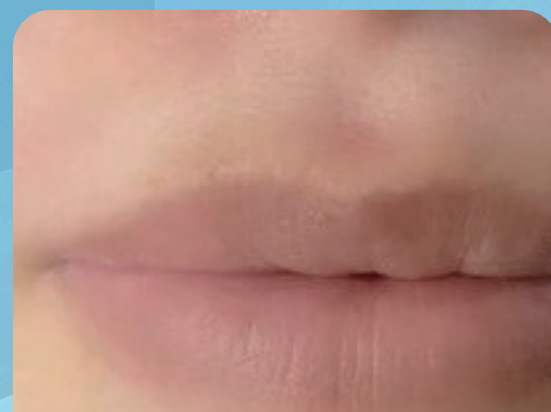

After

## Frequently Asked Questions (FAQ)

### HOW MANY SESSIONS ARE REQUIRED FOR HAIR REMOVAL?

Most patients see significant results after 6–8 sessions, spaced 4–6 weeks apart

### IS TRILUX ICE SAFE FOR DARKER SKIN TONES?

Yes, the device is safe for Fitzpatrick I–VI, thanks to customisable wavelengths and advanced cooling

### DOES THE TREATMENT HURT?

No, ICE cooling technology ensures treatments are pain-free, even for sensitive skin.

### HOW LONG DOES A SESSION TAKE?

- Small areas (e.g., upper lip): 5 minutes.
- Large areas (e.g., back or legs): 15–20 minutes.

# Device Specifications

- **Wavelength Options:** 755nm, 808nm, 1064nm.
- **Spot Size:** 12mm x 24mm for large areas
- **Frequency:** 1–10Hz for fast, efficient treatments
- **Pulse Width:** Adjustable 10–400ms
- **Cooling System:** TEC + water + air cooling for pain-free sessions
- **Power Output:** 600W–1200W for consistent results
- **Durability:** Tested for over 50 million shots

## Transform Your Clinic with TriLux ICE

Enhance your clinic's services with TriLux ICE, the ultimate diode laser for pain-free hair removal and skin rejuvenation. Its advanced technology and reliable performance make it the perfect investment for modern clinics.

### CONTACT US

✉ [sales@aestheticsconsults.com.au](mailto:sales@aestheticsconsults.com.au)

☎ +61 489 904212; +61 2 8526 4838

📍 Shop 23A Mt Annan Market Place 11-13 Main St, Mt Annan, NSW 2567



Contact Aesthetics Consults today to schedule a demo or learn more about TriLux ICE

